

VETPUMP® 2

Precision Suction &
Irrigation System



Cleaning ears has never been so easy

The VETPUMP® 2 from KARL STORZ puts the safest and most effective flushing and suction at your fingertips!

More powerful, more precise and safer than other units, the VETPUMP® 2 is the ultimate device for lavage and suction during otoscopy, gastroscopy and more.



Features:

- Powerful suction and irrigation pump for both rigid and flexible KARL STORZ endoscopes
- There is no exchange of air between bottles, eliminating any potential contamination of clean irrigation fluids
- Micro-processor controlled unit is quiet and responsive; unit is completely silent when trumpet valves are not activated
- Irrigation pressure and suction vacuum can be regulated separately
- User-selected pressures are precise, powerful and reproducible
- Versatile hand piece attaches to catheters, irrigation tubes, or directly to endoscope channel
- Bottles, tubing and handle are fully autoclavable

Check out our schedule of upcoming hands-on training courses at <http://go.karlstorz.com/eventsVET>





69 3216 01

69 3216 01

VETPUMP® 2, power supply 100/120/230/240 VAC, 50/60 Hz including:



Suction Bottle, 1.5 l, sterilizable
Bottle Cap, for use with suction bottle 1.5 l, sterilizable



Irrigation Bottle, 1 l, sterilizable
Bottle Cap, for use with irrigation bottle

Bottle Holder, for suction bottle
Bottle Holder, for irrigation bottle
Bacterial Filter

Mains Cord
Mains Cord, US-version

| Pressure Regulated (mmHg) | Flow Regulated (l/min) | Suction Pressure (mmHg) | Power Supply | Dimensions (w x h x d) mm | Weight (kg) | Certified to |
|---------------------------|------------------------|-------------------------|-------------------------------|---------------------------|-------------|--------------|
| max. 600 | max. 3.5 | regulated: max. 580 | 100/120/230/240 VAC, 50/60 Hz | 305 x 110 x 260 | 3.8 | EN 61010-1 |

Accessories for Otoscopy:



69320033 **Silicone Tubing Set**, irrigation-suction from unit to bottles, **autoclavable**



69811 **Instrument Tubing Set**, autoclavable, to connect VETPUMP® 2 bottles to handle



69810 **Handle**, for irrigation and suction, autoclavable, for use with tubing set 69811, in conjunction with 5 mm suction and irrigation tubes or catheters any size with LUER Lock connection



Accessories for Abdominal/Thoracic Lavage:



69320033 **Silicone Tubing Set**, irrigation-suction from unit to bottles, **autoclavable**



69811 **Instrument Tubing Set**, autoclavable, to connect VETPUMP® 2 bottles to handle



37112 A **Handle**, straight, with clamping valve, for suction and irrigation, **autoclavable**, for use with 5 and 10 mm suction and irrigation tubes



69811 LAP **Adaptor**, for use with suction and irrigation Handle 37112 A, plus Instrument Tubing Set 69811, package of 5



031121-10* **Gas filter**, for single use, sterile, package of 10, with male LUER-Lock connector for attachment to irrigation connector on unit



Accessories for Gastrointestinal Endoscopy:



20300042 **Silicone Tubing Set**, for suction only, sterilizable



610 VE **Silicone Tubing Set**, for insufflation, attaches to light source adaptors 487 VE or 69487 KVE




*




VETPUMP® 2

System Components for use in Otoscopy

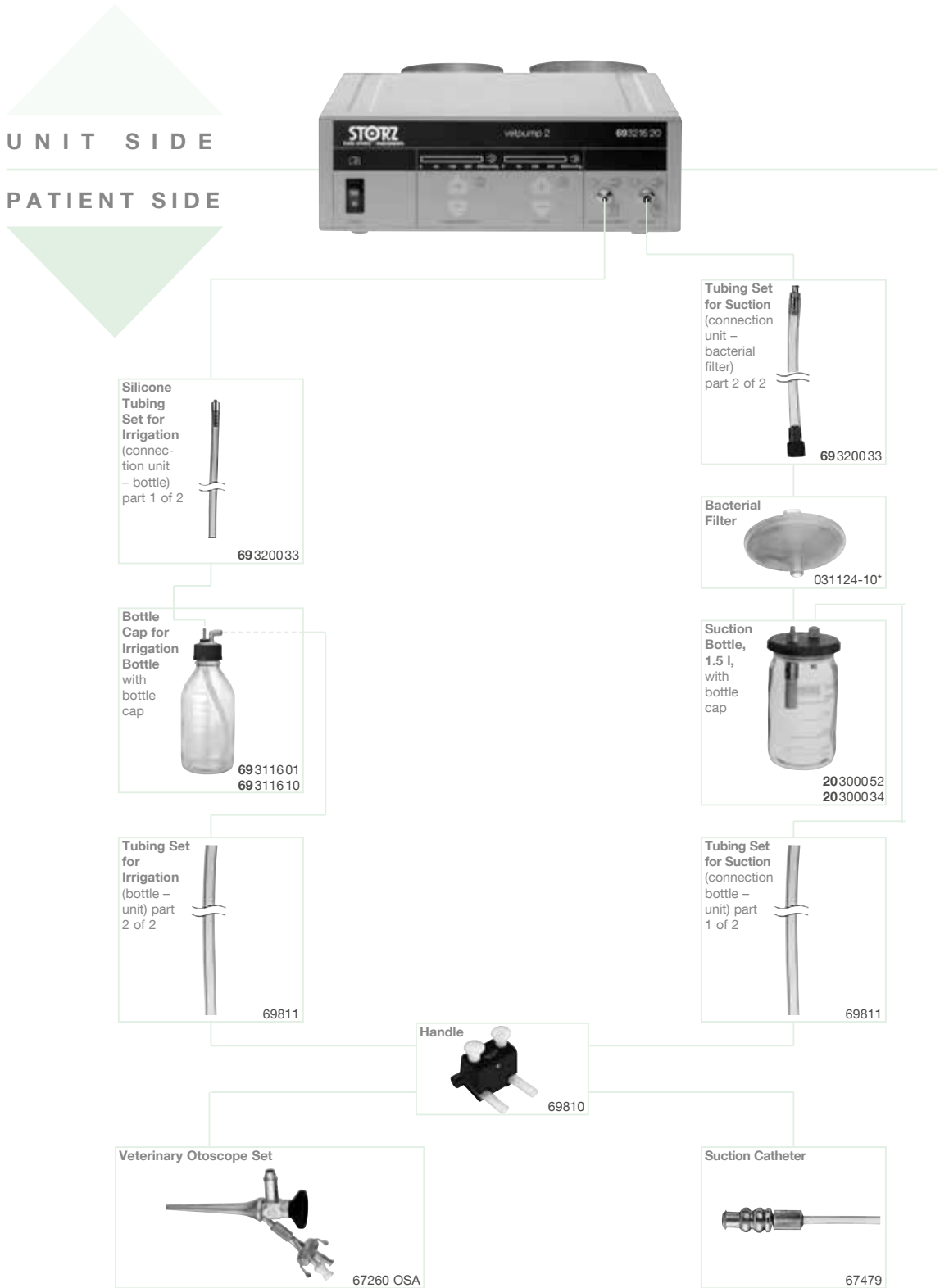

VETPUMP® 2

System Components for use in Laparoscopy / Thoracoscopy

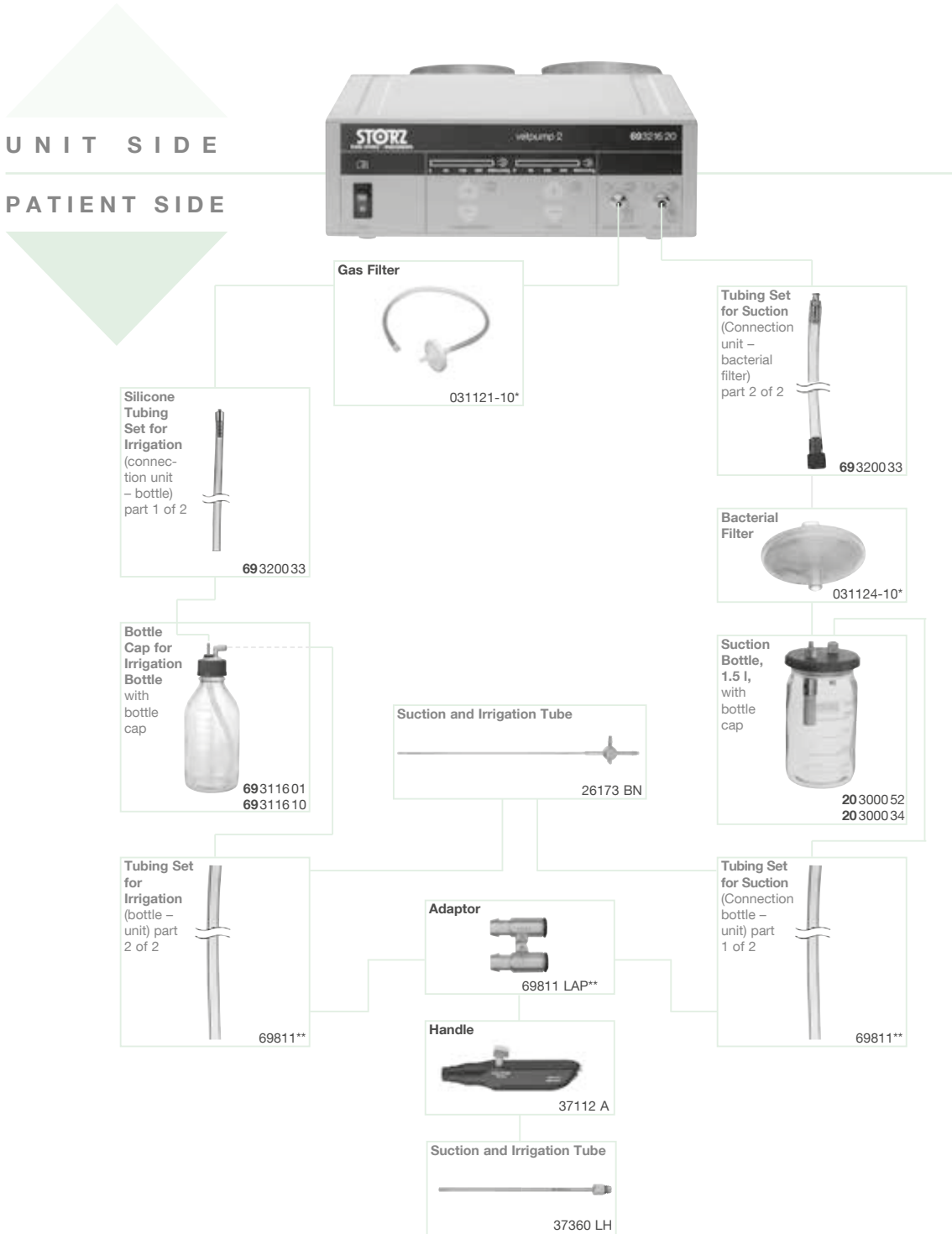


**Remove black plastic fitting from Tubing Set 69811 and replace with Adaptor 69811 LAP.

*



VETPUMP® 2

System Components for use in Gastrointestinal Endoscopy

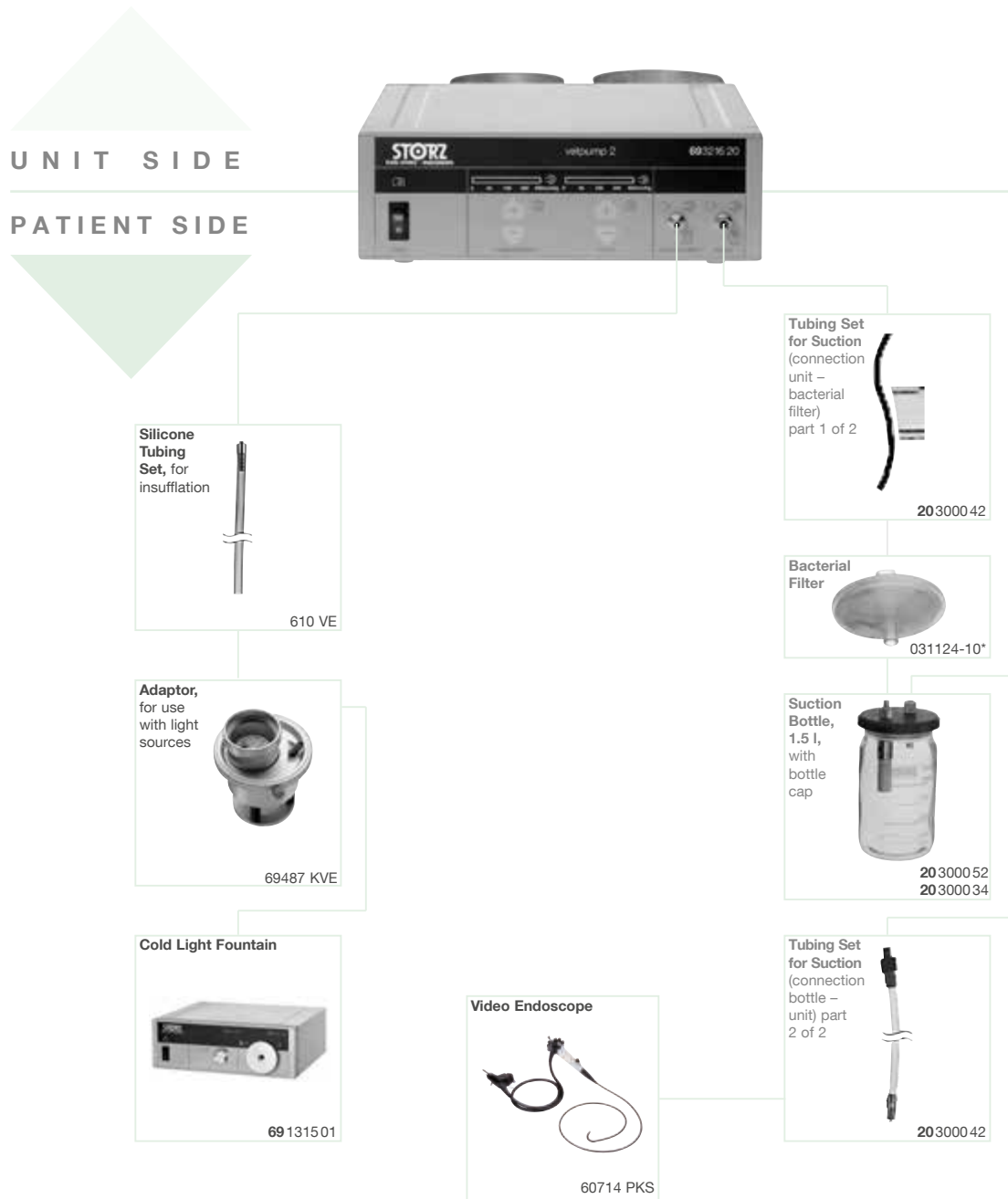

** Remove male LUER Lock fitting from end of Tubing **20300042** and connect free end to suction connector of endoscope.



KARL STORZ SE & Co. KG
Dr.-Karl-Storz-Straße 34, 78532 Tuttlingen/Germany
Postbox 230, 78503 Tuttlingen/Germany
Phone: +49 (0)7461 708-0
Fax: +49 (0)7461 708-105
E-Mail: vet@karlstorz.com
www.karlstorz.com