


C-MAC® PLATINUM SERIES AIRWAY MANAGEMENT CART

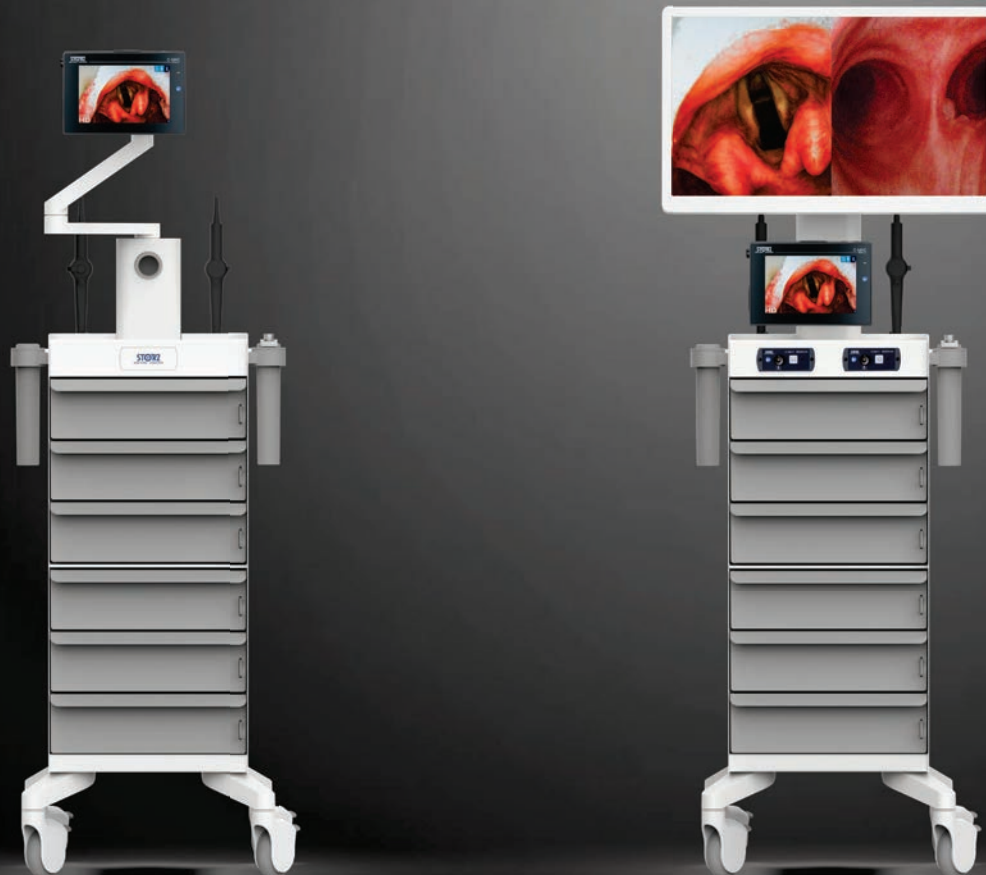

Platinum Cart Series

- Organization and secure storage of your video laryngoscopes, flexible scopes and accessories.
- C-MAC® devices close at hand; keeping you prepared for your “Plan B”.
- Featuring a lightweight design with adjustable monitor mounts for individual users.

C-MAC® HD
CHOICE CHANGES EVERYTHING

www.karlstorz.com

STORZ
KARL STORZ – ENDOSKOPE

C-MAC® PLATINUM SERIES AIRWAY MANAGEMENT CART

Choice Changes Everything!

The Platinum Carts offer multiple options to choose from for various storage needs. The Platinum Carts provide quick access to both video laryngoscopes and flexible intubation scopes, saving precious time when switching between airway devices. The smooth rolling lightweight Platinum Carts, provide storage for a comprehensive airway management platform at a great value.

- **Monitor Mount:** The Platinum Cart has a VESA mounting plate for attaching the C-MAC® monitor. The Platinum Plus has a stationary VESA mount that supports up to 19.8 lbs and a VESA mount for the C-MAC® monitor
- **Two Side Holders:** Provide storage for the S-Imager or Video Stylet, which can be easily removed for cleaning
- **Six Drawers:** Lined with the same Santoprene™ material, slide smoothly, and can be locked with red cable tie tags
- **Two Scope Holders:** Provide protection and storage for flexible intubation scopes. Includes a holster on the back to hold two boxes of FIVE S Single-use scopes
- **Articulating Arm:** The articulating arm on the Platinum Cart allows the C-MAC® monitor to tilt up/down and pan left/right, for maximum adjustability for individual users
- **Casters:** The precision-bearing casters provide a smooth and quiet experience, with locking casters for stability when static.



Side Holders



Two Scope Holders



Brackets for C-HUB II

Part Number: XUSM0011
Description: PlatinumPlus Cart for Video Intubation Platform
Dimensions: 18.9" (H) × 20.6" (W) × 59.9" (D)
Weight: 127 lbs.

Santoprene is a registered trademark of Exxon Mobil Corp

Part Number: XUSM0012
Description: Platinum Cart for Video Intubation Platform
Dimensions: 18.9" (H) × 20.6" (W) × 59.8" (D)
Weight: 105 lbs.

Santoprene is a registered trademark of Exxon Mobil Corp

**Contact your local KARL STORZ representative for more information
about these or other KARL STORZ products, or call 800-421-0837.**

It is recommended to check the suitability of the product for the intended procedure prior to use.
Please note that the described products in this medium may not be available yet in all countries due to different regulatory requirements.

KARL STORZ SE & Co. KG
Dr.-Karl-Storz-Straße 34, 78532 Tuttlingen/Germany
Postbox 230, 78503 Tuttlingen/Germany
Phone: +49 7461 708-0
Fax: +49 7461 708-105
E-Mail: info@karlstorz.com
www.karlstorz.com



KARL STORZ Endoscopy-America, Inc.
2151 East Grand Avenue
El Segundo, CA 90245-5017, USA
Phone: +1 424 218-8100
Phone toll free: 800 421-0837 (US only)
Fax: +1 424 218-8525
Fax toll free: 800 321-1304 (US only)
E-Mail: communications@karlstorz.com