




MTP Single-use Products for irrigation and suction in laparoscopy and gynecology



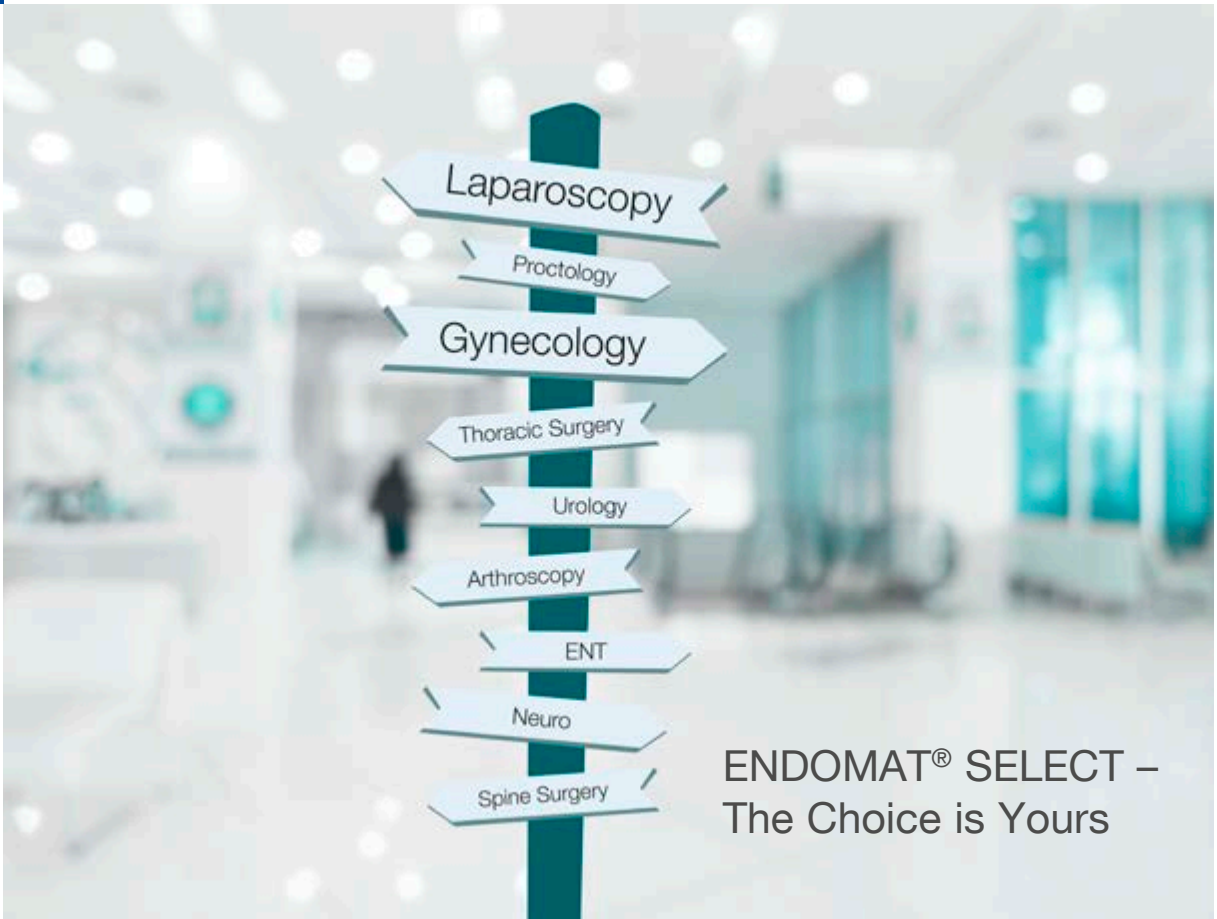

ENDOMAT® SELECT – The Choice is Yours

No matter what you choose – laparoscopy, gynecology, proctology, thoracic surgery, urology, arthroscopy, ENT, NEURO or spine surgery – the ENDOMAT® SELECT offers reliable support!


ENDOMAT® SELECT is a single-roller pump that is used for the irrigation or suction of fluids during surgical and diagnostic procedures. The pump can be equipped with various software modules so that the unit can be used for multiple disciplines in various rooms according to the respective indication.

Single-use Tubing Sets for use with ENDOMAT® SELECT


KARL STORZ offers indication-specific, single-use tubing sets for the ENDOMAT® SELECT. By clearly displaying the speciality on the packaging, the user is able to determine at a glance which tubing set is suitable for which indication. Thanks to a [coding](#) system, tubing cassettes with one-hand control are only recognized by the unit if the appropriate software has been installed. This means that the pump only functions within the parameters permitted for a specific application.

Tubing Cassettes Coding

Irrigation

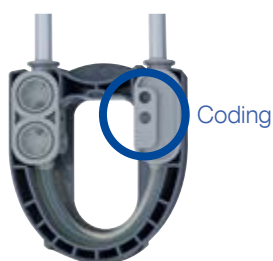


Flow control

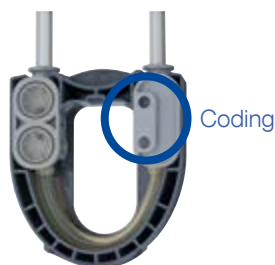

Pressure control

Suction

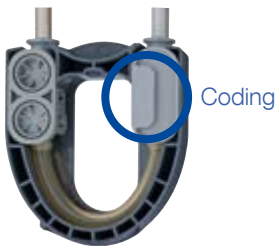

Direct suction

Benefits of using single-use vs. reusable tubing sets

In comparison to reusable tubing sets, the single-use tubing set avoids the risk of cross-contamination as the tubing set is only used once on a single patient.

This eliminates the need for time-consuming dismantling, pre-cleaning, drying, assembly, packaging and resterilization, thus minimizing the use of resources for the user. KARL STORZ single-use tubing sets are delivered sterile and are ready for immediate use.

Benefits of using a pump vs. working with gradient pressure

When working with gradient pressure, it is not possible to control the pressure or the flow. The use of the ENDOMAT® SELECT, on the other hand, offers the advantage of flow and/or pressure regulation – depending on the application.

Suction/Irrigation Handles for Single Use

For use with pump

The suction/irrigation handle for use with the pump is delivered with a suction tube but without an irrigation tube.

The irrigation tubing set* that is compatible with a particular pump is connected to the **tube olive of the handle**. Patient-side silicone tubes with an inner diameter of 5 mm can also be connected to the tube olive.

030138-10 **Irrigation/Suction Handle with Suction/Irrigation Probe****, sterile, for single use, package of 10



Tube olive

For use without pump

The suction/irrigation handle for use without a pump is suitable for work with gradient pressure and is delivered with both a suction and an irrigation tube.

030238-10 **Irrigation/Suction Handle with Suction/Irrigation Probe****, sterile, for single use, package of 10



*The irrigation tubing set is not included in delivery

**Diameter 5 mm, length 330 mm

Suction/Irrigation Handles for Multiple Use

KARL STORZ offers various models of reusable suction/irrigation handles. These handles are used in combination with Tubing Set 031134-10, which is equipped with a tube connector. The tubing set* compatible with the used pump is attached to this tube connector. Tubing Set 031133-10 is required for reusable suction/irrigation handles used without a pump.



37113A **Handle**, pistol grip, with clamping valve, for suction and irrigation, autoclavable



37112A **Handle**, straight, with clamping valve, for suction and irrigation, autoclavable

For use with pump



031134-10 **Tubing Set**, sterile, for single use, package of 10
STERILE 

*The irrigation tubing set is not included in delivery



37113A



37112A



031134-10

Tube connector



The irrigation tubing set that is compatible with the used pump and patient-side silicone tubes with an inner diameter of 5 mm can be attached to the [tubing connector](#).

For use without pump



031133-10

Tubing Set, sterile, for single use, package of 10



Overview of Single-use Tubing Sets for ENDOMAT® SELECT for use in Laparoscopy and Gynecology



031524-10 **Tubing Set**, for irrigation, FC, sterile, for single use, package of 10, for use with KARL STORZ HAMOU® ENDOMAT® 26331120-1 and ENDOMAT® SELECT UP210

STERILE 




031134-10 **Tubing Set**, for irrigation, sterile, for single use, package of 10, for use with KARL STORZ Handles 37112A (straight), 37113A (pistol grip) in combination with patient-side silicone tube (inner diameter 5 mm)


STERILE 



031133-10 **Tubing Set**, sterile, for single use, package of 10, for use with KARL STORZ Handles 37112A (straight), 37113A (pistol grip) in combination with patient-side silicone tube (inner diameter 5 mm)


STERILE 

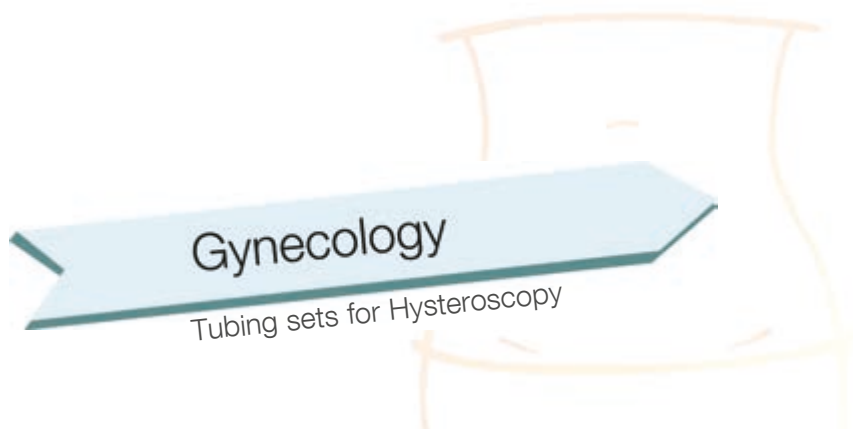

031523-10 **Tubing Set**, for irrigation, PC, sterile, for single use, package of 10, for use with KARL STORZ HAMOU® ENDOMAT® 26331120-1 and ENDOMAT® SELECT UP210
STERILE 





For combination with UNIDRIVE® S III

030647-10 **Tubing Set**, for suction, DS, sterile, for single use, package of 10, for use with KARL STORZ ENDOMAT® SELECT UP210
STERILE 



Now it's your choice!

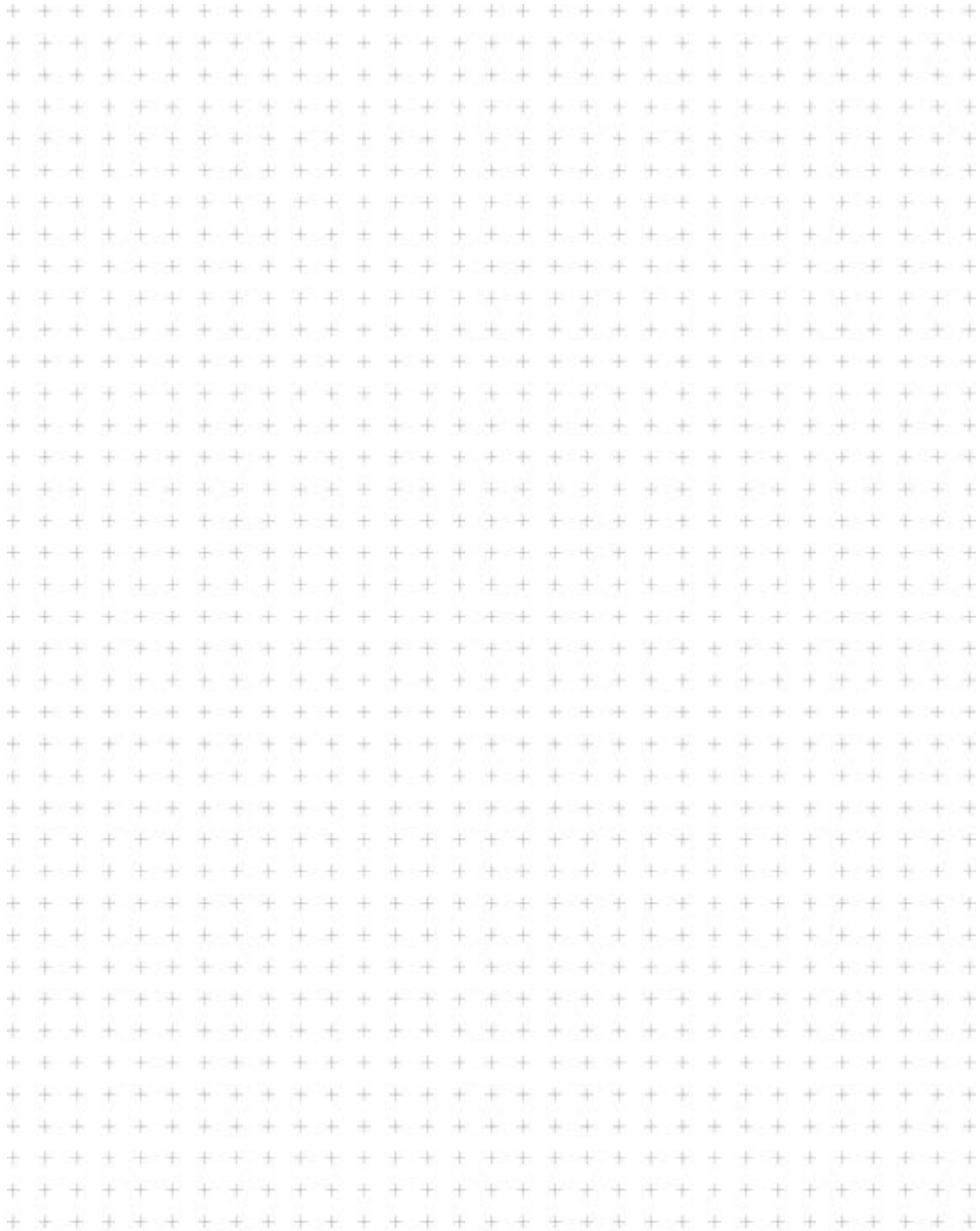

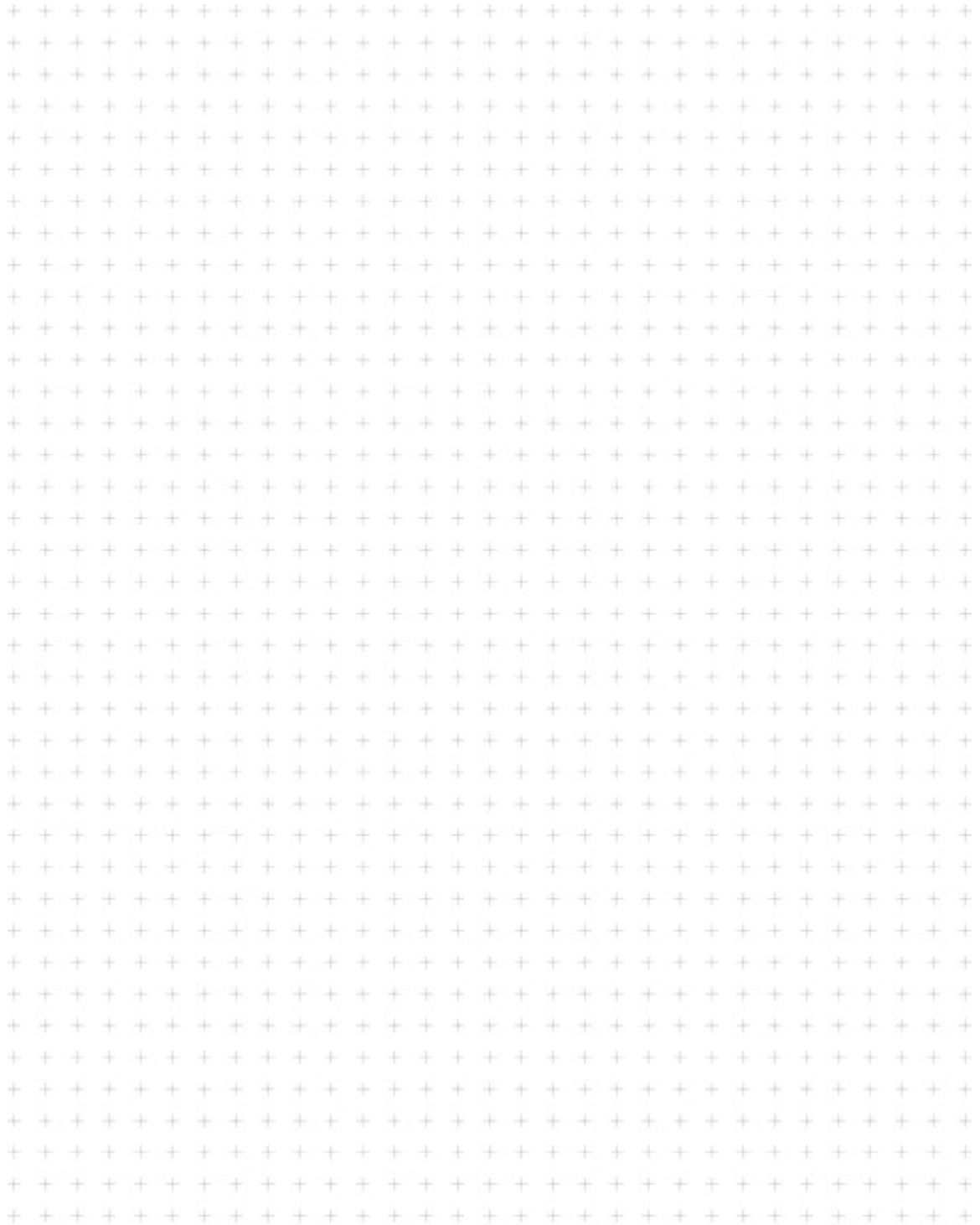

It is recommended to check the suitability of the product for the intended procedure prior to use.
Please note that the described products in this medium may not be available yet in all countries due to different regulatory requirements.

More than
75
Years

*Shaping the Future
of Endoscopy with you*

STORZ
KARL STORZ—ENDOSKOPE
THE DIAMOND STANDARD

KARL STORZ SE & Co. KG
Dr.-Karl-Storz-Straße 34, 78532 Tuttlingen/Germany
Postbox 230, 78503 Tuttlingen/Germany
Phone: +49 7461 708-0
Fax: +49 7461 708-105
E-Mail: info@karlstorz.com
www.karlstorz.com

